


visocolor® ECO Chlorid


Chloride / Chlorure / Cloruro

Messbereich / Range / Domaine de mesure / Rango
1 - 60 mg/l Cl⁻


REF 931 018
931 218

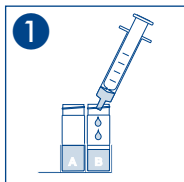


1
2 x 5 ml Probe
Sample
Echantillon
Muestra

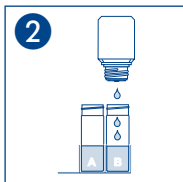


2
10 H_2O Cl-1

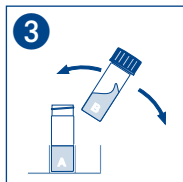


3
Schütteln
Shake
Agiter
Agitar

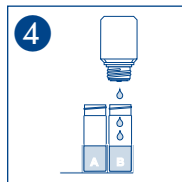


4
10 H_2O Cl-2

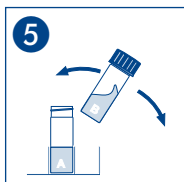

5
Schütteln
Shake
Agiter
Agitar

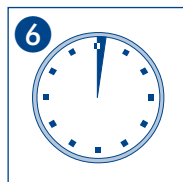

6
1'00 min

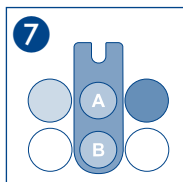

7
Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar