


# visocolor® ECO Cyanid


## Cyanide / Cyanure / Cianuro

Messbereich / Range / Domaine de mesure / Rango  
0.01 - 0.20 mg/l CN<sup>-</sup>


REF 931 022  
931 222

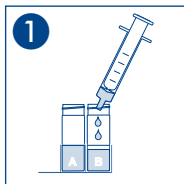

2 x 5 ml Probe  
Sample  
Echantillon  
Muestra

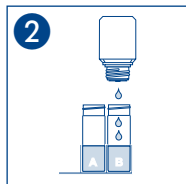


5 CN-1

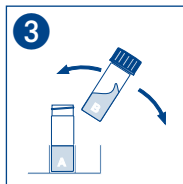


Schütteln  
Shake  
Agiter  
Agitar




1 CN-2

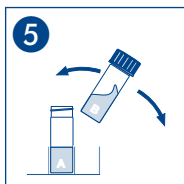


Umschwenken  
Invert  
Brandir  
Ondear

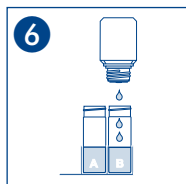


5 CN-3

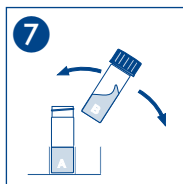


Schütteln  
Shake  
Agiter  
Agitar

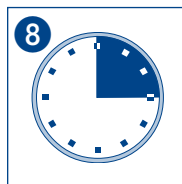

15'00 min

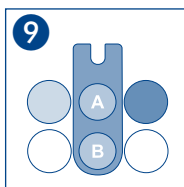

Messung  
Measurement  
Mesure  
Medición

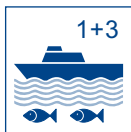

1+3

Meerwasser  
Sea water  
Eau de mer  
Agua de mar