


# visicolor® ECO Gesamthärte

## total Hardness / Dureté totale / Durezza total


Messbereich / Range / Domaine de mesure / Rango

1 °d  $\triangleq$  1 °d


REF 931 029

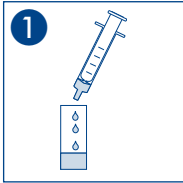


**1**  
5 ml Probe  
Sample  
Echantillon  
Muestra

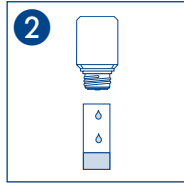


**2** ° GH-1

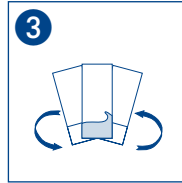

**3**  
Umschwenken  
Shake gently  
Secouer légèrement  
Mezclar volteándolo

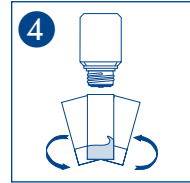

**4** ° GH-2  
bis Farbumschlag  
till coloration  
jusque inversion de couleur  
hasta cambio de color

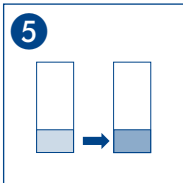

**5**  
1 °d  $\triangleq$  1 °d  
rot → grün  
red → green  
rouge → vert  
rojo → verde



1:30

Meerwasser  
Sea water  
Eau de mer  
Agua de mar