

visocolor® ECO Kieselsäure


Silica / Silice

Messbereich / Range / Domaine de mesure / Rango


0.2 - 3.0 mg/l SiO₂

REF 931 033


931 233

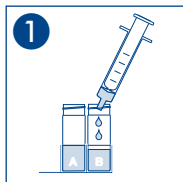


1
2 x 5 ml Probe
Sample
Echantillon
Muestra

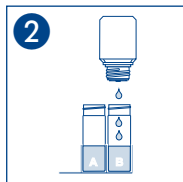


2
5 ⚡ SiO₂-1

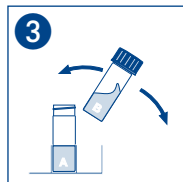


3
Schütteln
Shake
Agiter
Agitar

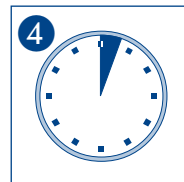


4
3'00 min

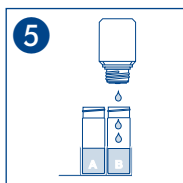


5
5 ⚡ SiO₂-2

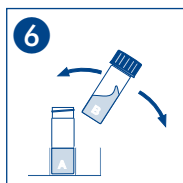


6
Schütteln
Shake
Agiter
Agitar

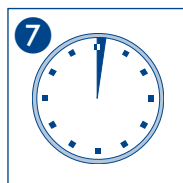


7
1'00 min

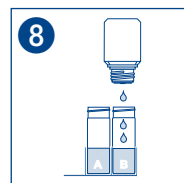


8
5 ⚡ SiO₂-3

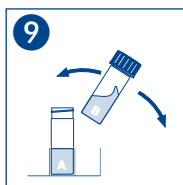

9
Schütteln
Shake
Agiter
Agitar

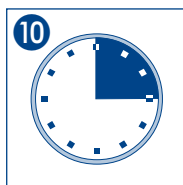

10
15'00 min

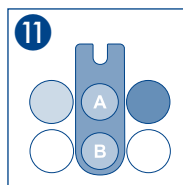

11
Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar