


visicolor® ECO Sulfat


Sulphate / Sulfate / Sulfato

Messbereich / Range / Domaine de mesure / Rango
25 - 200 mg/l SO₄²⁻


REF 931 092
931 292
914 035
914 235

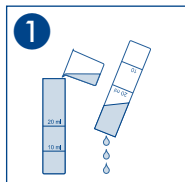


Spülen
Rinse
Rincer
Lavar

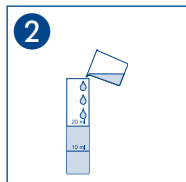


20 ml Probe
Sample
Echantillon
Muestra

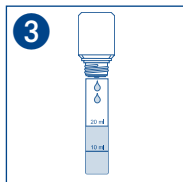


10 Δ SO₄⁻¹

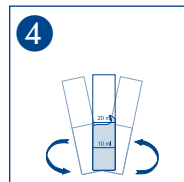


Umschwenken
Shake gently
Secouer légèrement
Mezclar volteándolo

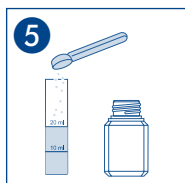


1 Δ SO₄⁻²

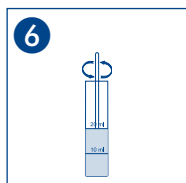


Umrühren
Stir
Remuer
Revolver

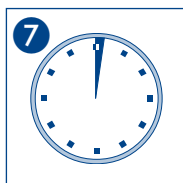

1'00 min

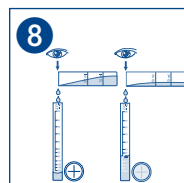

Messröhrchen füllen
Filling measuring tube
explétif la cuve de mesure
echar en tubito de medida

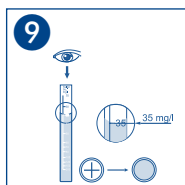

bis Kreuz unsichtbar
till cross invisible
jusque croix invisible
hasta cruz invisible



1:50

Meerwasser
Sea water
Eau de mer
Agua de mar