


visocolor® ECO Sulfid

Sulphide / Sulfure / Sulfuro



Messbereich / Range / Domaine de mesure / Rango
0.1 - 0.8 mg/l S²⁻


REF 931 094
931 294

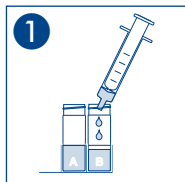


1
2 x 5 ml Probe
Sample
Echantillon
Muestra




2
1 S-1

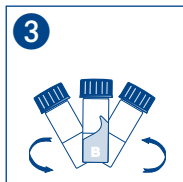


3
Umschwenke
Shake gently
Secouer légèrement
Mezclar volteándolo

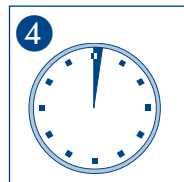


4
1'00 min

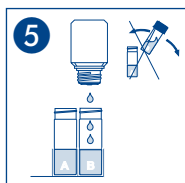


5 S-2

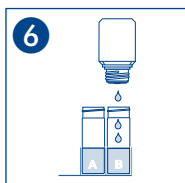


6 S-3

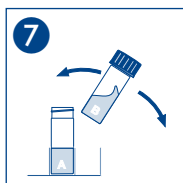

7
Schütteln
Shake
Agiter
Agitar

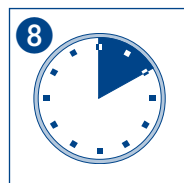

8
10'00 min

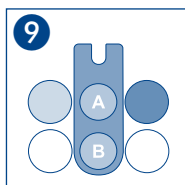

9
Messung
Measurement
Mesure
Medición



Meerwasser
Sea water
Eau de mer
Agua de mar