


visicolor® ECO Zink

Zinc




Messbereich / Range / Domaine de mesure / Rango
0.5 - 3 mg/l Zn²⁺


REF 931 098
931 298

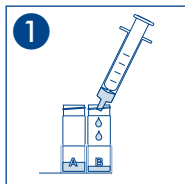


1
2 x 1 ml Probe
Sample
Echantillon
Muestra

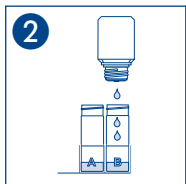


2
2 Zn-1

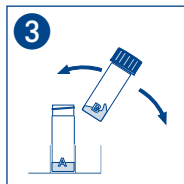


3
Schütteln
Shake
Agiter
Agitar

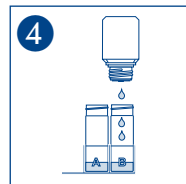


4
2 Zn-2

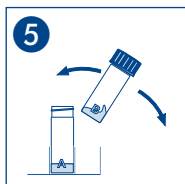


5
Schütteln
Shake
Agiter
Agitar

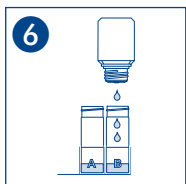


6
5 Zn-3

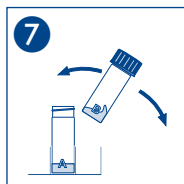


7
Schütteln
Shake
Agiter
Agitar

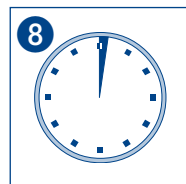

8
1'00 min

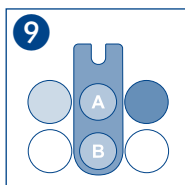

9
Messung
Measurement
Mesure
Medición

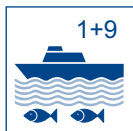

1+9

Meerwasser
Sea water
Eau de mer
Agua de mar