




Fluoride / Fluorure / Fluoruro

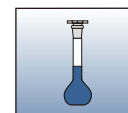


Cat. No.	Cuvette	Measurement range	Method	Wavelength
918 142	10 mm	0.05 - 2.00 mg/l F ⁻	(1) 421	585 nm
	50 mm	-		

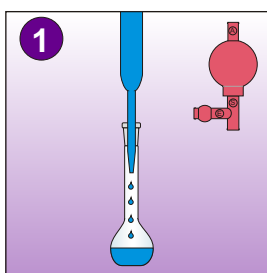


1 20 ml sample
 échantillon
 muestra

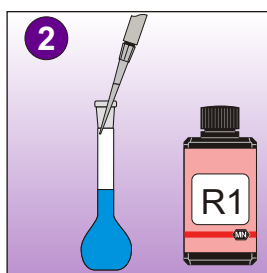


2 2.0 ml R1

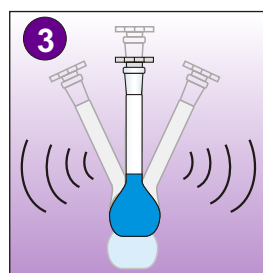


3 mix
 mélanger
 mezclar

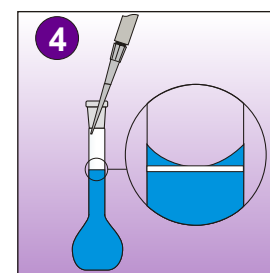


4 fill up with dist. water
 emplir avec l'eau dist.
 llenar con agua dest.

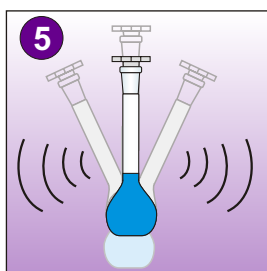


5 mix
 mélanger
 mezclar

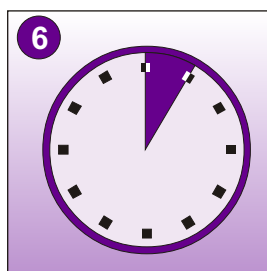


6 5'00 min

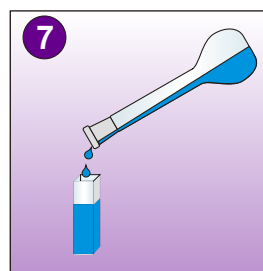


7 rinse and fill
 rincer et emplir
 lavar y llenar

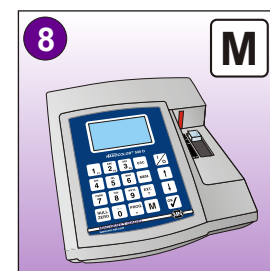

8 measurement
 mesure
 medición

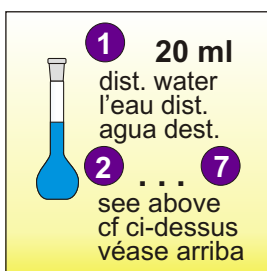

blank value
 zéro
 cero

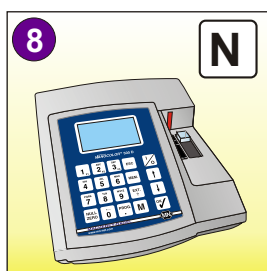

blank measurement
 zéro mesure
 cero medición