




# Chlorine (free/total) / Chlore / Cloro

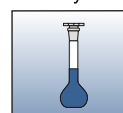


Cat. No.	Cuvette	Measurement range	Method	Wavelength
918 16	10 mm	0.1 - 10.0 mg/l Cl <sub>2</sub>	free (1) 161 total (1) 162	540 nm
	50 mm	0.02 - 2.00 mg/l Cl <sub>2</sub>	free (1) 161 total (1) 162	

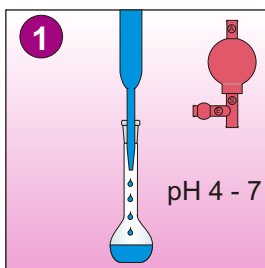


**1** 20 ml sample  
échantillon  
muestra

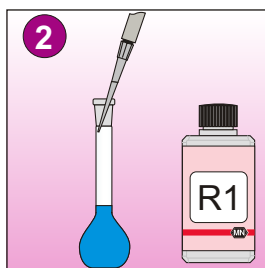


**2** 1.0 ml R1

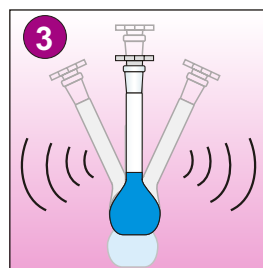


**3** mix  
mélanger  
mezclar

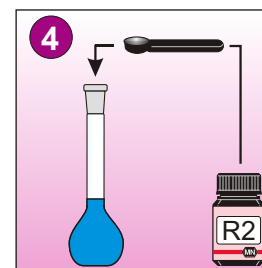

**4** 1 R2

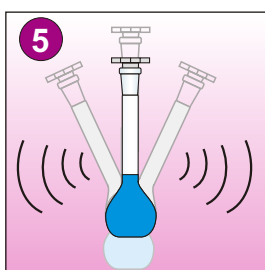


**5** mix  
mélanger  
mezclar

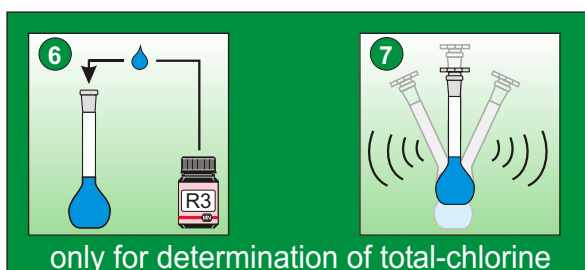


**6** 5 R3


**7** mix  
mélanger  
mezclar

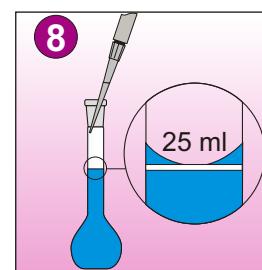

**8** fill up with dist. water  
emplir avec l'eau dist.  
llenar con agua dest.

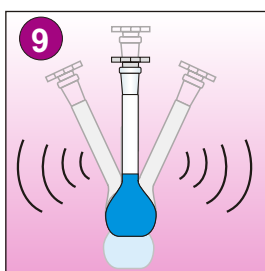


**9** mix  
mélanger  
mezclar

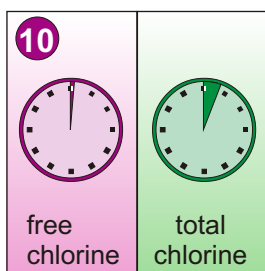


**10** 30 s


**10** 3 min

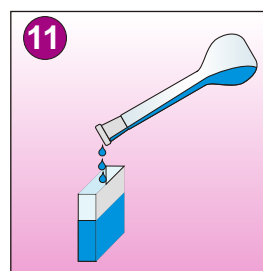


**11** rinse and fill  
rincer et emplir  
lavar y llenar

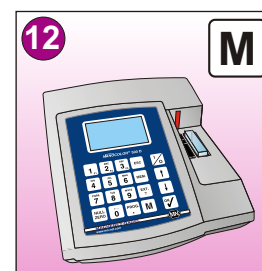

**12** measurement  
mesure  
medición

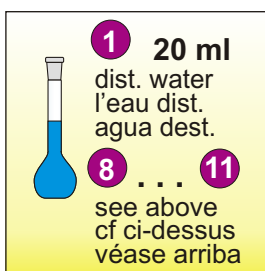

blank value  
zéro  
cero

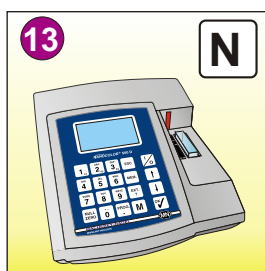

blank measurement  
zéro mesure  
cero medición