




# Cyanide / Cyanure

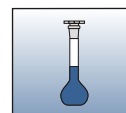


Cat. No.	Cuvette	Measurement range	Method	Wavelength
918 30	10 mm	0.01 - 0.50 mg/l CN <sup>-</sup>	(1) 301	585 nm
	50 mm	0.001 - 0.100 mg/l CN <sup>-</sup>		

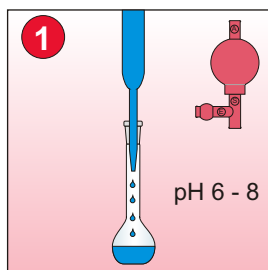


**1**  
20 ml sample  
échantillon  
muestra

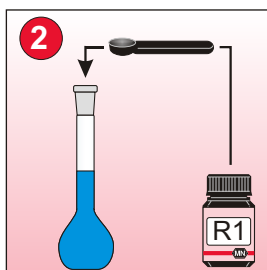


**1** R1

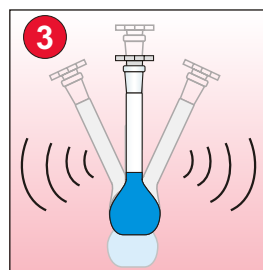


**3**  
mix  
mélanger  
mezclar

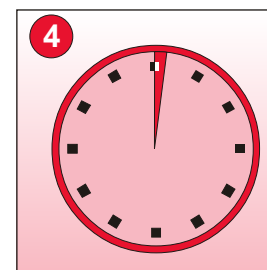


**4**  
1'00 min

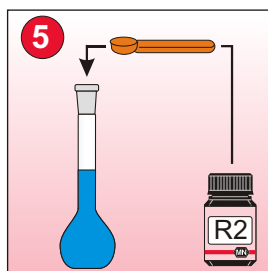


**1** R2

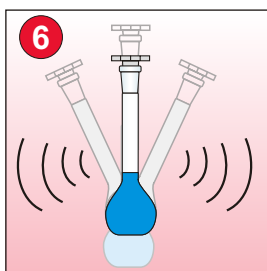


**6**  
mix  
mélanger  
mezclar

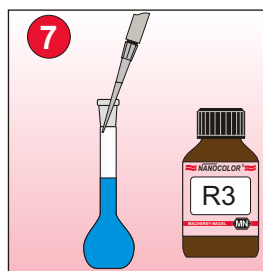


**7**  
2.0 ml R3

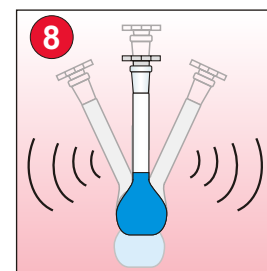


**8**  
mix  
mélanger  
mezclar

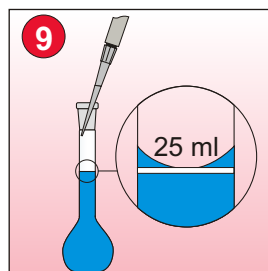


**9**  
fill up with dist. water  
emplir avec l'eau dist.  
llenar con agua dest.

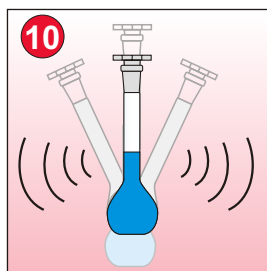

**10**  
mix  
mélanger  
mezclar

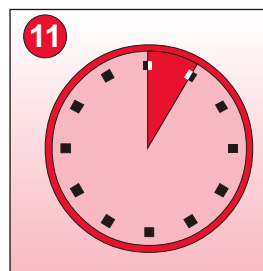


**11**  
5'00 min

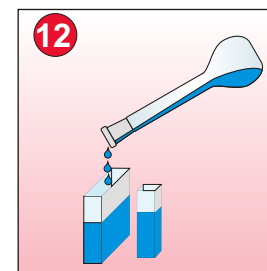

**12**  
rinse and fill  
rincer et emplir  
lavar y llenar



**13** M  
measurement  
mesure  
medición

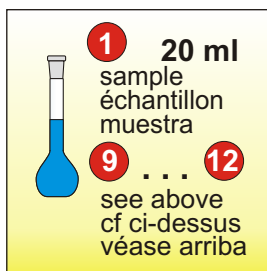

**13** N  
blank value  
zéro  
cero



**13** N  
blank measurement  
zéro mesure  
cero medición