




Detergents cation / Détergents cation

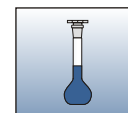


Cat. No.	Cuvette	Measurement range	Method	Wavelength
918 34	10 mm	0.2 - 5.0 mg/l CTAB	(1) 341	436 nm
	50 mm	0.05 - 3.00 mg/l CTAB		

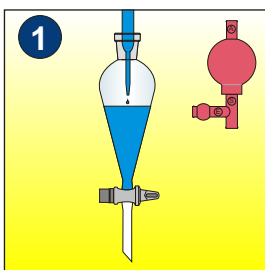


1
50 ml sample
échantillon
muestra

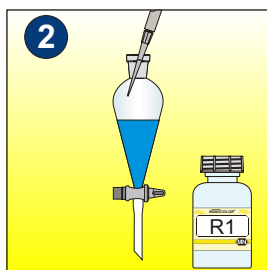


2
2.0 ml R1

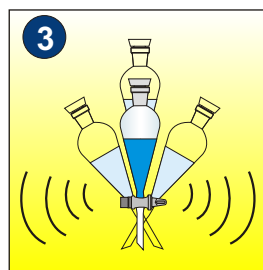


3
mix
mélanger
mezclar

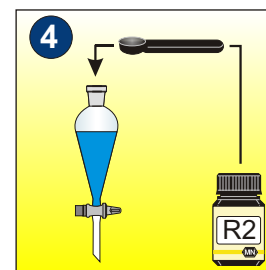


4
1 ml R2

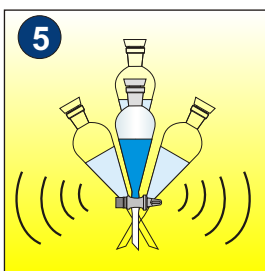


5
mix
mélanger
mezclar

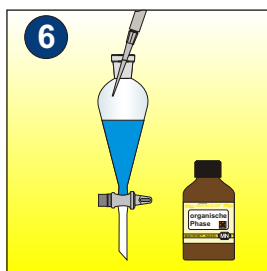


6
20 ml org. phase
org. fase

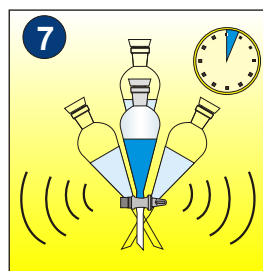

7
3 min shake
agiter
agitar

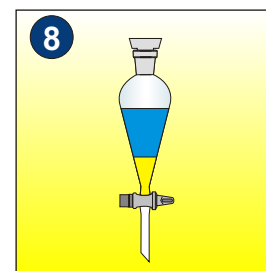


8
allow to separate
séparer
separar

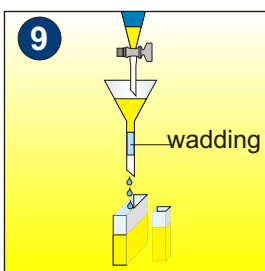


9
filtrate into cuvette
filtrer en cuves
filtrar en cubeta



10 **M**
measurement
mesure
medición

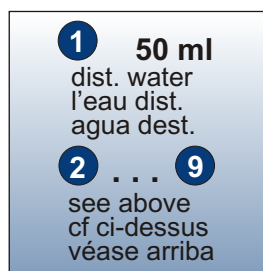

10 **N**
blank value
zéro
cero

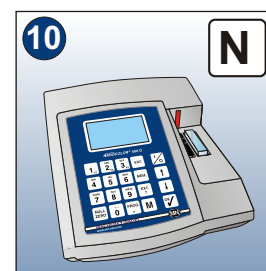

10 **N**
blank measurement
zéro mesure
cero medición