




Iron / Fer / Hierro

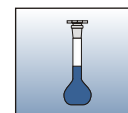


Cat. No.	Cuvette	Measurement range	Method	Wavelength
918 36	10 mm	0.1 - 15.0 mg/l Fe	(1) 361	470 nm
	50 mm	0.01 - 2.00 mg/l Fe		

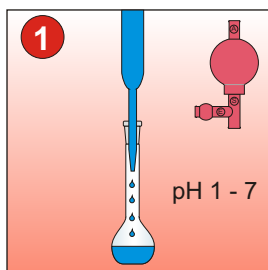


1
 20 ml sample
 échantillon
 muestra

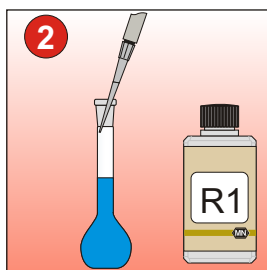


2
 1.0 ml R1

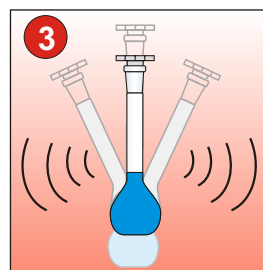


3
 mix
 mélanger
 mezclar

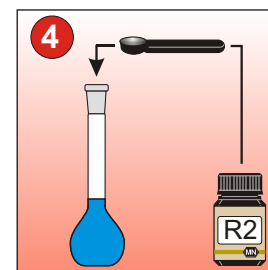


4
 1 spoon R2

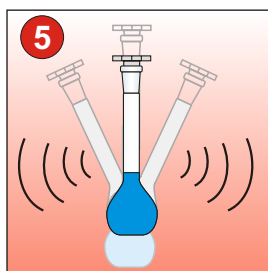


5
 mix
 mélanger
 mezclar

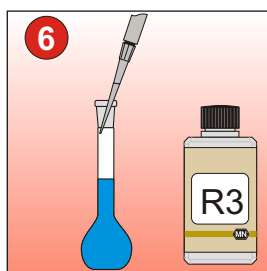


6
 1.0 ml R3

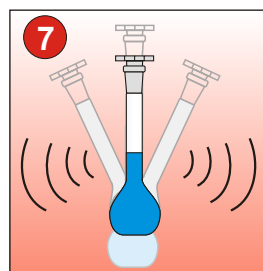


7
 mix
 mélanger
 mezclar

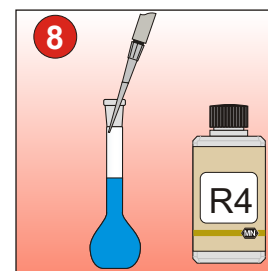


8
 1.0 ml R4

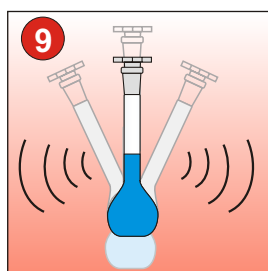


9
 mix
 mélanger
 mezclar

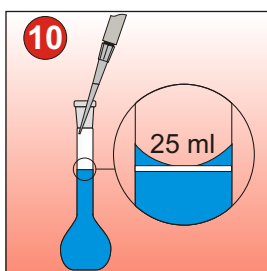


10
 fill up with dist. water
 emplir avec l'eau dist.
 llenar con agua dest.

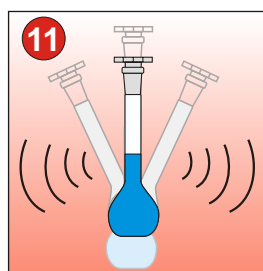

11
 mix
 mélanger
 mezclar

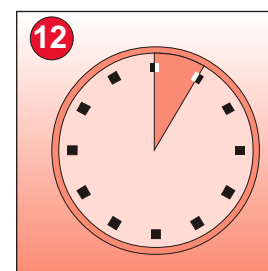


12
 5'00 min

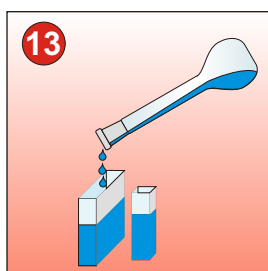


13
 rinse and fill
 rincer et emplir
 lavar y llenar



14 **M**
 measurement
 mesure
 medición

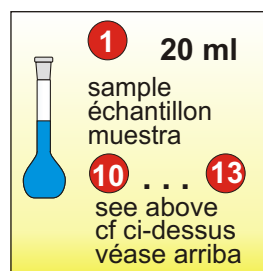

14 **N**
 blank value
 zéro
 cero

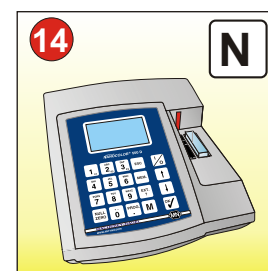

14 **N**
 blank measurement
 zéro mesure
 cero medición

1 20 ml
 sample
 échantillon
 muestra
10 ... **13**
 see above
 cf ci-dessus
 véase arriba