




Manganese / Manganèse / Manganeseo

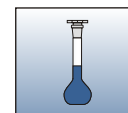


Cat. No.	Cuvette	Measurement range	Method	Wavelength
918 60	10 mm	0.1 - 10.0 mg/l Mn	(1) 601	470 nm
	50 mm	0.01 - 2.00 mg/l Mn		

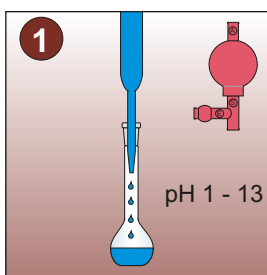


1
20 ml sample
échantillon
muestra

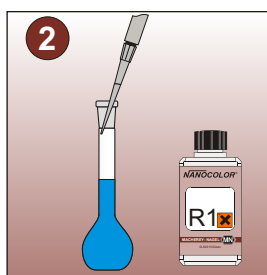


2
1.0 ml R1

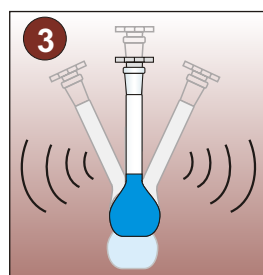


3
mix
mélanger
mezclar

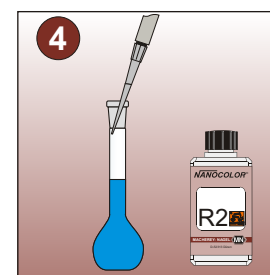


4
1.0 ml R2

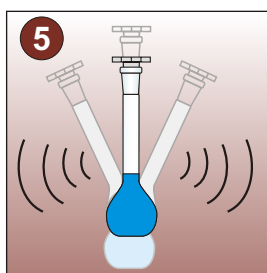


5
mix
mélanger
mezclar

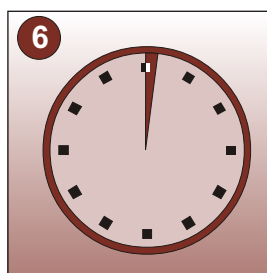


6
1'00 min

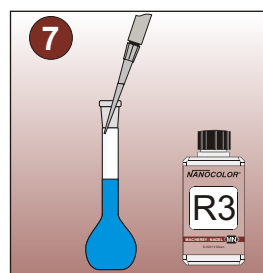


7
1.0 ml R3

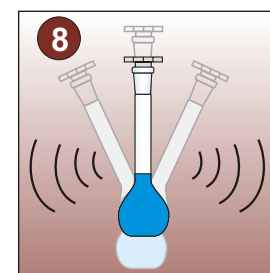


8
mix
mélanger
mezclar

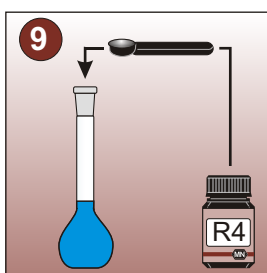


9
1 ml R4

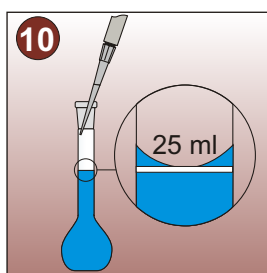


10
fill up with dist. water
emplir avec l'eau dist.
llenar con agua dest.

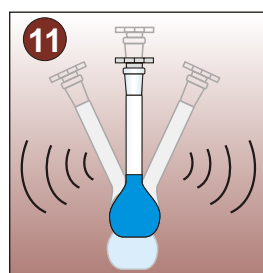


11
mix
mélanger
mezclar

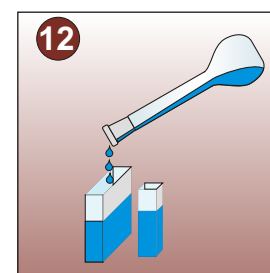


12
rinse and fill
rincer et emplir
lavar y llenar

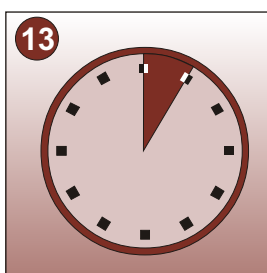


13
5'00 min

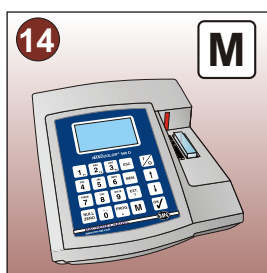

14
measurement
mesure
medición

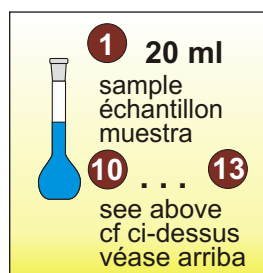

14
blank value
zéro
cero

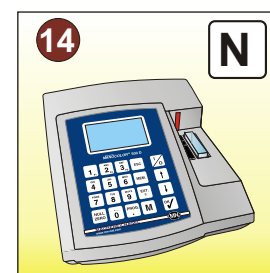

14
blank measurement
zéro mesure
cero medición