




Nitrite / Nitrito

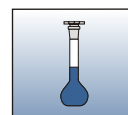


Cat. No.	Cuvette	Measurement range	Method	Wavelength
918 67	10 mm	0.01 - 0.30 mg/l NO ₂ ⁻ - N 0.10 - 1.00 mg/l NO ₂ ⁻	(1) 671 (1) 672	520 nm
	50 mm	0.002 - 0.100 mg/l NO ₂ ⁻ - N 0.005 - 0.250 mg/l NO ₂ ⁻	(1) 671 (1) 672	

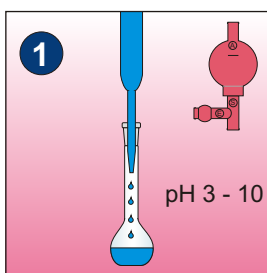


1
20 ml sample
échantillon
muestra

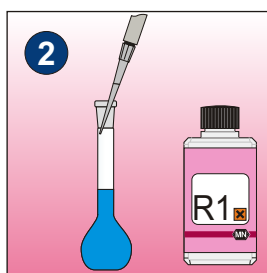


2
2.0 ml R1

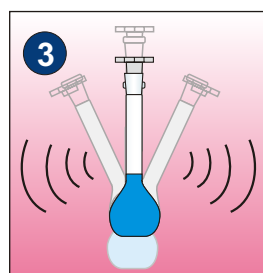


3
mix
mélanger
mezclar

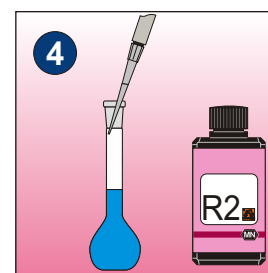


4
2.0 ml R2

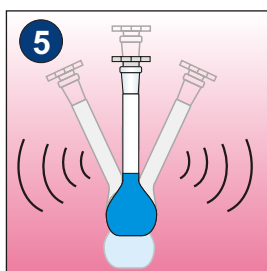


5
mix
mélanger
mezclar

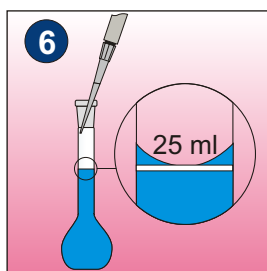


6
fill up with dist. water
emplir avec l'eau dist.
llenar con agua dest.

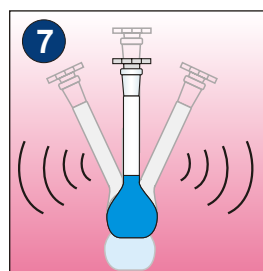


7
mix
mélanger
mezclar

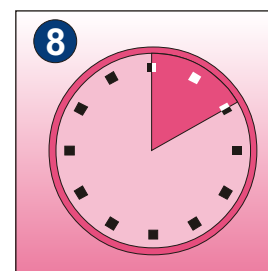


8
10'00 min

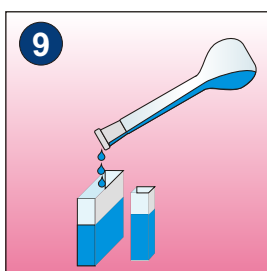


9
rinse and fill
rincer et emplir
lavar y llenar

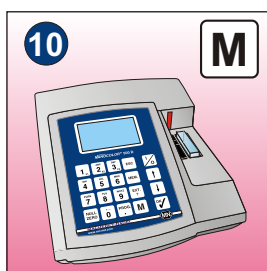

10
measurement
mesure
medición

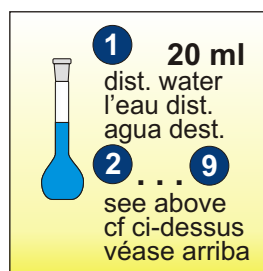

10
blank value
zéro
cero

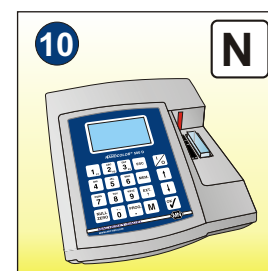

10
blank measurement
zéro mesure
cero medición