




Phenol / Phénol / Fenol

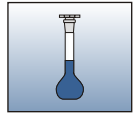


Cat. No.	Cuvette	Measurement range	Method	Wavelength
918 75	10 mm	0.1 - 7.0 mg/l Phenol	(1) 751	470 nm
	50 mm	0.01 - 1.00 mg/l Phenol		

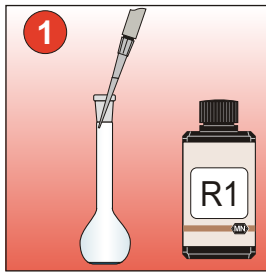


1.0 ml R1

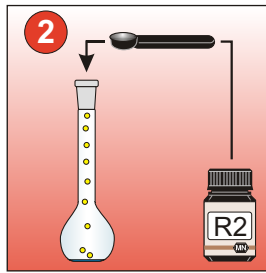


a few crystals R2 til decolorization
 quelque particules R2 jusqu'à décoloration
 unas cuantas particulas R2 hasta decoloración

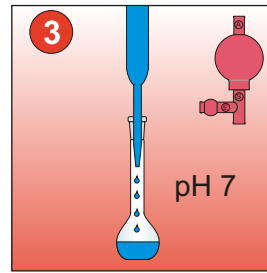


20 ml sample
 échantillon
 muestra

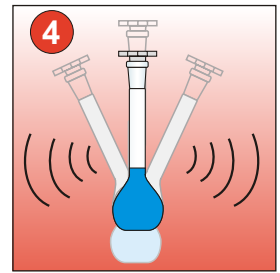


mix
 mélanger
 mezclar

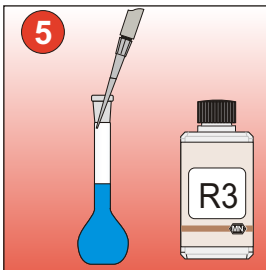


1.0 ml R3

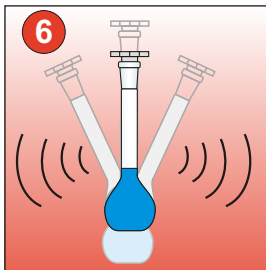


mix
 mélanger
 mezclar

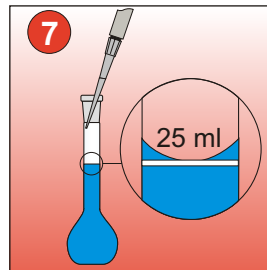


fill up with dist. water
 emplir avec l'eau dist.
 llenar con agua dest.

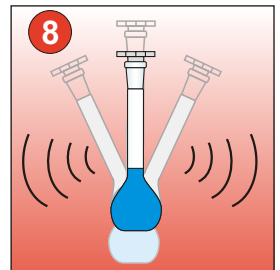

mix
 mélanger
 mezclar

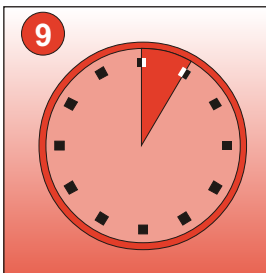


5'00 min

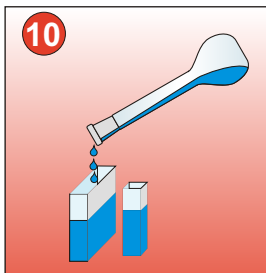


rinse and fill
 rincer et emplir
 lavar y llenar



measurement
 mesure
 medición

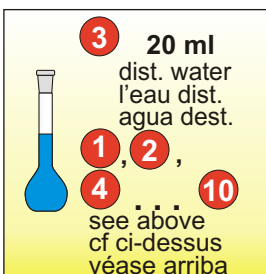

blank value
 zéro
 cero

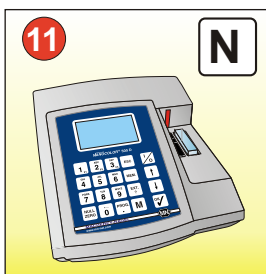

blank measurement
 zéro mesure
 cero medición