




ortho-Phosphate / ortho-Fosfato

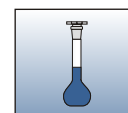


| Cat. No. | Cuvette | Measurement range | Method | Wavelength |
|----------|---------|--|-------------------------------|------------|
| 918 78 | 10 mm | 1 - 17 mg/l PO ₄ - P 2 - 50 mg/l PO ₄ ³⁻ 1.5 - 37.0 mg/l P ₂ O ₅ | (1) 781 (1) 782 (1) 783 | 436 nm |
| | 50 mm | 0.2 - 6.6 mg/l PO ₄ - P 0.5 - 20.0 mg/l PO ₄ ³⁻ 0.4 - 15.0 mg/l P ₂ O ₅ | (1) 781 (1) 782 (1) 783 | |

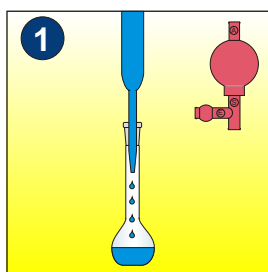


1
20 ml sample
échantillon
muestra

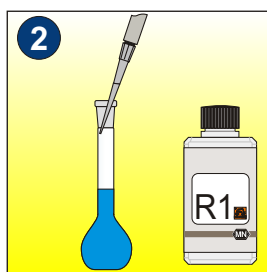


2
1.0 ml R1

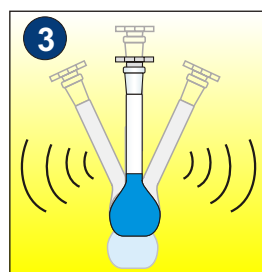


3
mix
mélanger
mezclar

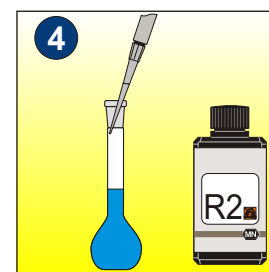


4
1.0 ml R2

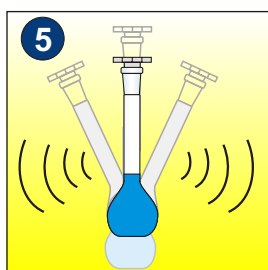


5
mix
mélanger
mezclar

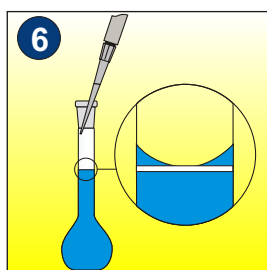

6
fill up with dist. water
emplir avec l'eau dist.
llenar con agua dest.

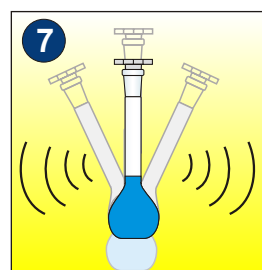

7
mix
mélanger
mezclar

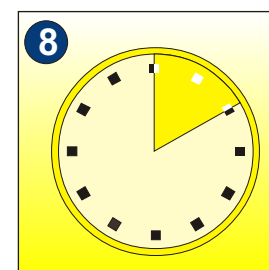


8
10'00 min




9
rinse and fill
rincer et emplir
lavar y llenar



10 **M**
measurement
mesure
medición

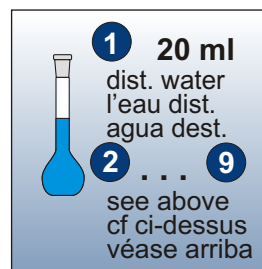

10 **N**
blank value
zéro
cero

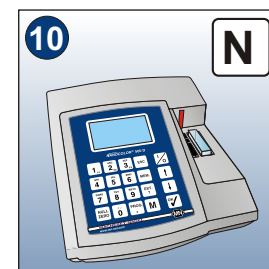

10 **N**
blank measurement
zéro mesure
cero medición